

Tourniquet System

Ideal for single or bi-lateral orthopaedic surgery and pain management



Innovative **Medical** Technology - **Practically Applied**


AneticAid

AT4™ TOURNIQUET SYSTEM

Intuitive and simple to use.

Compact, light weight and robust, offering advanced technology combined with maximum safety and reliability.



The AT4™ Electronic Tourniquet is a truly independent device operating from its own integral air supply, so no need for air supply hoses, medical air cylinders or power leads.

Engineered to use minimum power, the 25 hour battery can be in use all day without the need to recharge until the end of the list. The battery provides more than back-up: it's independence and reliability.

The Tourniquet has sensitive leak compensation facilities as well as automatic self-diagnostic checks on start-up, which means it is calibrated every time it is switched on. Cuff channel hoses and displays are all colour coded for ease of recognition.

SAFE - EFFECTIVE - PROVEN

Accessories

As well as the tourniquet machine itself, we offer a full range of accessories including three types of tourniquet cuff:

- REUSABLE** traditional, reusable and robust cuffs for long-term use
- SIX-USE** an ideal alternative to re-usable cuffs offering reduced risk of cross-infection
- SINGLE-USE** providing sterility when the risk of infection is key

These are complemented by foam cuff liners and cuff covers which protect the patient's limbs from pinching and shear forces when the tourniquet is inflated.



Your AT4 Tourniquet comes with an optional 5 year manufacturer's warranty.

This is subject to our usual Terms and Conditions around regular servicing which can be obtained through AnetiCare – planned, routine maintenance programmes which offer safety and peace of mind through a range of simple, structured pricing plans.

Call 01943 878647 to find out more.



The AT4 Electronic Tourniquet:

A state-of-the-art, dual channel electronic tourniquet with integral power supply, which requires no external compressed air supply and comes with mobile stand and cuff supply hoses. Anetic Aid® are the distributor of the AT4™ Electronic Tourniquet which is manufactured by Medical Device Management Ltd., 31 Braintree Business Park, Blackwell Drive, Braintree, CM7 2PU.

Description	Code
AT4 Electronic Tourniquet System	40080



TOURNIQUET FEATURES

Designed to be superbly manoeuvrable, stable and reliable



Automatic self-diagnostic checks
On start-up and throughout use.

IVRA (Bier's block)
Procedural guidance feature

Easy to clean
Low maintenance

Automatic LCD timers
With user defined elapsed time alarm

Dual channel system
Single or dual channel operation

Colour coded dual channel system
From controls to cuff - single or dual channel operation

Quick release hosing
For ease of storage

Truly independent
Device operates from its own integral air supply, so no need for air supply hoses, medical air cylinders or power leads

Bright OLED Independent displays
Pre-select and actual cuff pressure readings

Visual and audible warnings
For battery, leaks and calibration



Precision controlled regulators
User-friendly, simple to use controls

Additional storage facility

Positive feel push buttons
Inflate / deflate functions

Pressure range
0 – 600 mmHg

Cuff hooks
For ease of storage



Internal power supply
Rechargeable battery that lasts all day – not just battery backup



Combining industry knowledge with innovative technology



Designed and produced at our UK manufacturing facility using the latest engineering techniques



Manufacturer, designer and supplier of high quality medical equipment around the world

SIX-USE DAY TOURNIQUET CUFFS

Six-use Day Tourniquet Cuffs are an ideal alternative to re-usable cuffs, offering reduced risk of cross infection.

Guaranteed for six uses, they are cost effective enough to dispose of as a single use item. Latex free with a slim, low profile design, they are available in a range of colour coded sizes and supplied in boxes of 10.

Full Range of Sizes

Colour coded for easy recognition

Cost Effective

Less time on cleaning / checking for re-use

Light Weight

Also latex free

Slim, Low Profile

Designed for patient comforts

Secure Fastening

Eliminates rolls or twists



Six-use Day Tourniquet Cuffs (box of 10 – Straight)

Description	Code
■ Six-use Day Tourniquet Cuff – Paediatric	40285
■ Six-use Day Tourniquet Cuff – Adult Arm Narrow	40288
■ Six-use Day Tourniquet Cuff – Adult Arm	40290
■ Six-use Day Tourniquet Cuff – Small Leg Narrow	40292
■ Six-use Day Tourniquet Cuff – Small Leg	40293
■ Six-use Day Tourniquet Cuff – Adult Leg	40295
■ Six-use Day Tourniquet Cuff – XL Adult Leg	40296

Six-use Day Tourniquet Cuffs (box of 10 – Conical)

Description	Code
■ Six-use Day Tourniquet Cuff – Conical Adult Leg	40297
■ Six-use Day Tourniquet – XL Conical Adult Leg	40298

Patients come in all shapes and sizes and in some cases conical cuffs conform more closely to the contours of the limb, giving a better tourniquet block.

Exsanguinators

Description	Code
Rhys Davies Exsanguinator	40400
Esmarch Surgical Bandage 5.5m – Latex Free	40450

Rhys Davies Exsanguinator

The Rhys Davies Exsanguinator is a simple but effective piece of equipment which allows the exsanguination of limbs for surgery without subjecting the limb to the high pressures experienced with other methods of exsanguination.



Tourniquet Cuff Covers (box of 30)

Description	Code
Tourniquet Cuff Cover – Padded – Small	40305
Tourniquet Cuff Covers – Padded – Large	40315



TOURNIQUET CUFFS



Reusable Tourniquet Cuffs

Manufactured from easy-to-clean, robust, latex free materials with additional safety straps for long-term use, available in a variety of colour coded sizes from paediatric to extra-large adult leg.

Description	Code
Reusable Tourniquet Cuff – Paediatric	40200
Reusable Tourniquet Cuff – Arm	40220
Reusable Tourniquet Cuff – Arm – Narrow	40225
Reusable Tourniquet Cuff – Large Arm / Small Leg	40240
Reusable Tourniquet Cuff – Large Arm / Small Leg – Narrow	40245
Reusable Tourniquet Cuff – Large Arm / Small Leg	40240
Reusable Tourniquet Cuff – Leg	40260
Reusable Tourniquet Cuff – Leg – Conical	40268
Reusable Tourniquet Cuff – Large Leg	40270
Reusable Tourniquet Cuff – Large Leg – Conical	40278

Reusable Tourniquet Cuffs – IVRA / Bier’s Block

Description	Code
Reusable Tourniquet Cuff – Paediatric	40200
Reusable Tourniquet Cuff – Arm	40230
Reusable Tourniquet Cuff – Large Arm / Small Leg	40250

Anetic Aid offers a number of double tourniquet cuffs designed specifically for Bier’s block / IVRA procedures.

ADDITIONAL SAFETY	coated strap for long-term use
SLIM, LOW PROFILE	designed for patient comfort
SECURE FASTENING	eliminates roll or twist
LIGHT WEIGHT	also latex free
FULL RANGE OF SIZES	colour coded for ease of recognition



Single-use Sterile Tourniquet Cuffs (box of 10 - Straight)

Single-use, Sterile Tourniquet Cuffs (complete with stocking sleeve padding) are the latest addition to our range, providing sterility when the risk of infection is key.

Because they are sterile, these cuffs can be used in front of sterile drapes, giving the surgeon optimum operating space – particularly useful in paediatrics or for procedures such as total knee replacement or ulna nerve decompression.

Description	Code
Single-use Tourniquet Cuff – Sterile – Paediatric	40800
Single-use Tourniquet Cuff – Sterile – Adult Arm	40820
Single-use Tourniquet Cuff – Sterile – Small Leg	40840
Single-use Tourniquet Cuff – Sterile – Adult Leg	40860
Single-use Tourniquet Cuff – Sterile – XL Adult Leg	40870

Single-use Sterile Tourniquet Cuffs (box of 10 – Conical)

Description	Code
Single-use Tourniquet Cuff – Sterile – Conical Adult Leg	40868
Single-use Tourniquet Cuff – Sterile – XL Conical Adult Leg	40878

STERILE, SINGLE-USE	less time on cleaning / checking for re-use
SLIM, LOW PROFILE	designed for patient comfort
SECURE FASTENING	eliminates roll or twist
LIGHT WEIGHT	also latex free
FULL RANGE OF SIZES	colour coded for ease of recognition



WARRANTY
Subject to annual servicing



OUR SERVICE PLANS Protect your investment

We offer a unique and popular ‘5 Year Warranty’. Take out an **AnetiCare** service package for regular and predictable maintenance costs that includes a ‘5 Year Warranty’.



OUR SERVICE TEAM Nationwide coverage

Our team of manufacturer-trained technicians service your equipment on-site and will get to know your hospital or clinic inside out.



OUR EXPERTISE 5 year warranty

Covers engineer’s time and most parts. Contact us to learn more about what an **AnetiCare** service packaging can offer you. Available via your Anetic Aid partner - found at aneticaid.com/distributors



Value for money: minimum downtime, maximum peace of mind



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

SEE
ALSO



QA4 Mobile Surgery System:

The complete multi-purpose solution for transport, treatment and recovery

PLUS



our other product areas



Scan this code to start a WhatsApp chat with Anetic Aid



OUR PRODUCT AREAS

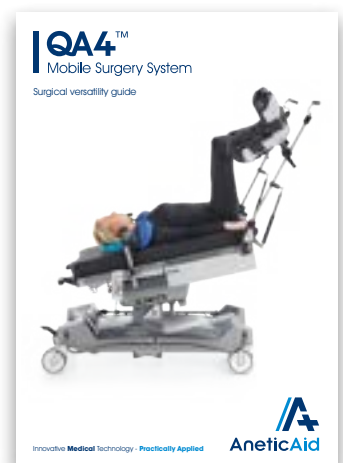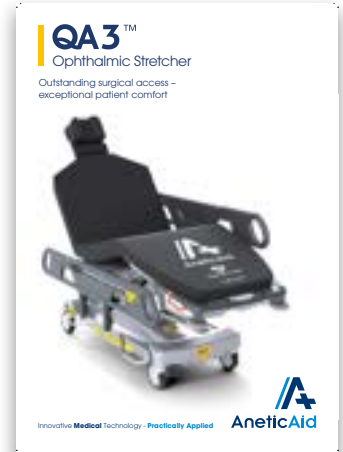

For more information on any of our products or service contracts, call:
sales@aneticaid.com +44 (0)1943 878647 aneticaid.com



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Innovative **Medical** Technology - **Practically Applied**



ISSUE 3 - FEB 2021