


QA3™

«800» Stretcher

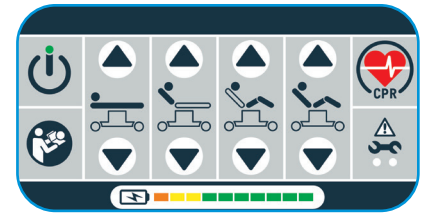

Powered ⚡ backrest

Positions the backrest anywhere between the vertical and horizontal for maximum patient comfort.

Inbuilt motorised //DRIVE™

No need to push or pull the stretcher: our new, inbuilt motorised DRIVE feature linked to the stretcher's 5th wheel, offers controlled, effortless travel up to gradients of 10°.



Patient comfort

Exceptional low height capability and store away side rails allow easy mount/dismount particularly for frail and elderly patients. Wide patient platform and pressure care mattress for patient comfort. Optimum ground clearance for patient hoist access.



Powered ⚡ platform height

Offers smooth movement to position the stretcher at any point in its height range - 475mm to 775mm (excluding mattress).

Castors

Polyurethane tyres offer low rolling resistance for high manoeuvrability, 150mm castors fitted as standard, 200mm optional.

Electrically powered functions

Easy to operate using the positive feel push button controls on a panel at the foot of the stretcher. Manual override systems as standard. Minimises risks for practitioners dealing with increasing patient weights.

SPECIFICATIONS			
Minimum height	*475 mm	Mattress length	2025 mm
Maximum height	*775 mm	Mattress width	800 mm
*excludes mattress		Mattress depth	105 mm
Minimum length	1700 mm	Maximum patient weight	320kg
Maximum length	2150 mm	Backrest	0 to +90°
Overall width	905 mm	KneeFlex	0 to +30°
Weight	128 kg	Trendelenburg	0 to -12°

sales@aneticaid.com +44 (0)1943 878647 aneticaid.com

Innovative **Medical** Technology - **Practically Applied**

Anetic Aid Ltd, Sapper Jordan Rossi Business Park, Otley Road, Baildon, BD17 7AX, UK



QA3™

◀800▶ Stretcher

The QA3 < 800 > is a special Bariatric version of the original QA3 with increased platform/mattress width.

◀800▶ for Bariatric Patients



All models of the QA3 Patient Stretcher have a maximum weight capacity of 320kg.



As with all of our stretchers, the QA3 < 800 > is fitted with the **K8 pressure care mattress**.



High density foam together with a layer of thermally responsive **gel-feel memory foam**.



Helps to **relieve pressure** by distributing the patient's weight evenly across the surface area.



Antibacterial and microbial fungistatic properties.



A **latex-free**, 4-way stretch polyurethane (PU) outer cover allows the skin to breathe - for **greater patient comfort**.

“ Our unique ultrasonic sealing system effectively welds the seams of the outer cover material together to give an almost unbreakable seal, making the K8 a class leader in infection control! ”